




FTM-Lite VIS • Film Thickness Gauge

High-Precision Film Thickness Desktop Gauge

The **FTM-Lite VIS** film thickness gauge uses the white-light interference phenomenon for non-contact and non-destructive film thickness determination of transparent single or double layers in the range of approximately 1 to 25 micrometer. The FTM-Lite VIS instrument is designed for manually performed, but even though high-precise measurements in the lab.

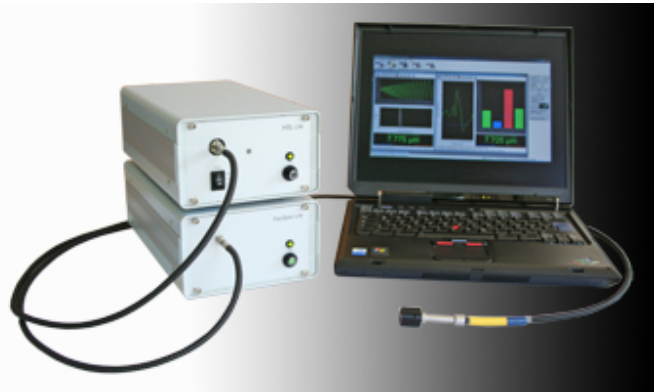

TranSpec Lite VIS - Desktop PDA-Spectrometer

Our TranSpec Lite VIS spectrometer exclusively uses the high-innovative spectrometer modules of Carl Zeiss, Germany. With these modules, the entrance slit is imaged on a photodiode array (PDA) by means of a holographic created, concave diffraction grating. All optical components of the spectrometer module are firmly mounted together in one unit and permanently adjusted to each other, which means that there are no mechanically moveable parts at all. Thus, the TranSpec Lite VIS spectrometer is free-of-maintenance!

HSL Lite - Desktop Halogen Spectral Lamp

The HSL Lite halogen spectral lamp is designed to operate with the TranSpec Lite VIS spectrometer in order to detect white-light interference spectra using a bifurcated fiber optics cable for your non-contact film thickness measurements. It comes with an integrated shutter, which can be manually opened and closed or automatically switched by TranSpec Lite VIS.

FTM-ProVis Lite - Powerful Film Thickness Software

This very easy-to-use software package uses an improved Fast-Fourier Transformation (FFT) algorithm to determine the film thickness from measured white-light interference spectra of thin transparent layers, which allows high-precision results in the entire measurement range. The film thickness result is computed in real-time, can be displayed in various different on-line charts and logged to a text file during the measurement. FTM-ProVis Lite also permits the fully automatic and simultaneous measurement of double-layers.

Technical specifications on next page ►



FTM-Lite VIS Film Thickness Gauge • Technical Specifications

May 2011, all information without guarantee and subject to changes

Minimum Hardware and Software Requirements

- PC/Laptop with Intel Dual-Core (2 GHz or higher recommended)
- Windows XP or Windows Vista or Windows 7 (32 bits version, native 64 bits in preparation)
- CD-ROM drive for installation
- Graphics adapter with at least 1024 x 768 pixel, 1280 x 1024 pixel is recommended
- One available USB 2.0 / 3.0 port

Mechanical Construction

- Desktop chassis with CE certificate, no cooling required
- Dimensions (H x W x D): - TranSpec Lite VIS: 100 x 160 x 280 mm - HSL Lite: 100 x 160 x 280 mm
- Weight - TranSpec Lite VIS: 1.8 kg - HSL Lite: 1.8 kg
- External power supply with 100/240 Volt auto-switch

TranSpec Lite VIS - Desktop PDA-Spectrometer

- Industry standard FSMA connector for bifurcated fiber optics cable
- Spectrometer modules of Carl Zeiss, Germany, with holographic created concave diffraction grating
- Photodiode array (PDA) with 256 pixel, no cooling required
- Permanently adjusted module, no mechanically moveable parts, maintenance free!
- Module specific wavelength range: ~ 300 - 1100 nm
- Module specific spectral resolution: ~ 10 nm
- Module specific spectral pixel interval: ~ 3.3 nm
- Absolute wavelength accuracy: ≤ 0.3 nm
- Temperature drift: typically < 0.005 nm / Kelvin

HSL Lite - Desktop Halogen Spectral Lamp

- Industry standard FSMA connector for bifurcated fiber optics cable
- Integrated mechanical shutter, manually or automatically switched
- 7 Watt miniature halogen bulb of Welch Allyn, Inc. USA
- Spectral output range: approximately 400 - 1000 nm
- Typical lifetime of the bulb: 1500 hours

FTM-ProVis Lite - Film Thickness Software

- Multi-threaded MDI-Application for Windows XP/Vista/7, online-help, color printed user's manual
- Evaluation of interference spectra with the help of a special Fast-Fourier-Transformation (FFT)
- New algorithm for high-precise sub-pixel determination of the FFT peak position (film thickness result)
- Accuracy typically within +/- 0.005 microns, repeatability typically within +/- 0.002 microns (standard WEG test)
- Real-time chart representation of Interference, FFT-Spectrum and Film Thickness Trend during measurement
- Consideration of refraction index and dispersion (Cauchy dispersion correction)
- Simultaneous measurement and film thickness evaluation of double layers possible

Film Thickness Measurement Range

The film thickness range generally measurable with FTM-Lite VIS is approximately 1 to 25 micrometer (~ 0.04 to 1 mil), but depends essentially on the currently selected spectral evaluation range, which can be setup in the software individually for each type of layer. Other factors which determine the measurable film thickness range are the refraction index (and its dispersion) of the layer to be measured.

Note TranSpec is a registered German trademark of Dipl.-Ing. (FH) Th. Fuchs, Engineer's Office for Applied Spectroscopy. All other mentioned product names are or possibly might be trademarks or registered trademarks of their owners.